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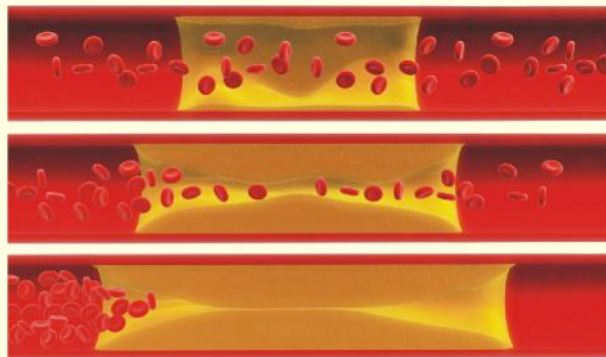
Operation Manual

【 Yu Yang Cup 】

According to the Master Handbook of Chinese Medicine “The Emperor’s Scripture” : Solar energy inside our body is the same as the energy of the Earth and the Sun. If you lose it, the longevity of life will be affected. Solar energy is vital, moves inside the body and strengthens the body. Yu Yang Cup is based on the theory of Chinese Medicine: Moxibustion (Heat) complements instinct energy and scrapping the skin improves the movements of meridians inside the body. Yu Yang Cup is a scientific product. It hands down the essence of Chinese Medicine, creates the nourishment and modernizes the way of keeping health.

While using Yu Yang Cup, different kinds of energy like moxibustion, warm scrapping, massage, magnetic treatment, heat treatment and infra-red energy can be simultaneously applied to the body, thus helps the solar energy to circulate inside the body. Good circulation not only improves the immune system to protect the body but also resists external harms and diseases.

【 Efficiency 】

Combined with the traditional Chinese medicine and modern medicine technology, the Yu Yang Cup can improve immune system, prevent from illness, and prolong life if regular use. It is indicated that by warming and scrapping of relevant acupuncture points can achieve the assistant treatment effects on relief of headache, insomnia, frozen shoulder pains, stiff neck pains, the waist and back pain. It can also help reduce the disorders of intestines and stomachs, dysmenorrheal and chronic exhaustion syndrome etc.

【 Functions of the product 】

1.Moxibustion: The product can produce heat varies from 45 degree Celsius to 85 degree Celsius to warm the meridians, unblock collaterals, disperse cold and eliminate wet . By means of the above procedures the operation of human organs can be in good conditions.

2.Warm scrapping: Scrapping the skin is also called “guasha” in Chinese Medicine. Scrapping the skin not only activates the collaterals but also removes blood stasis. Chinese Medicine Theory insists on “Good blood circulation can improve the immune system so as to prevent many kinds of illnesses.”

3.Infrared Ray: Infrared ray is also called “the light of life” , it helps improve blood circulation, softens the blood vessels, promotes nutritious status of tissues, relieves or eliminates pains, and improves immune system.

4.Magnetic therapy: According to scientific researches, magnetism can improve metabolism, activate cell regeneration effectively, eliminate inflammation and balance endocrine system.

【 How to operate 】

1.Turn on the main switch, connect the electric wire to the switch and plug in the Cup, choose the suitable temperature variable from 45°C to 85°C. The Cup will reach the preset temperature in about 5 minutes.

Air conditioning: The room temperature is better to be maintained at around 23°C with good ventilation.

2.By pressing the vibration button, the vibration of the Cup begins, a green light can be seen on the surface of the top.

3.It is easy to control the temperature of the Cup. Press “+” sign to increase the temperature to the most of 85°C and “-” to the least of 45°C. The default temperature is 65°C and every press varies 5°C.

4.The waveform at the bottom of the Cup is to scrape the muscles which are alongside with the bones and the governing meridian. The protrude part is used to apply on the trigger points.

5.The teeth form at the bottom is to scrape the stiff muscles at the back. It is a very good way to reduce the stress of the muscles.

【 Ways to protect Yu Yang Cup 】

Please use the cleaning fluid provided by our Company to clean the Cup after use. Just spray the fluid on the silica gel surface and the ceramic cap and clean the parts gently with dry tissue. Do not immerse it in water to avoid electric shock or any harm to the Cup. In case of any damage, do not use the Cup but immediately contact our company for repair. If the Cup is not to be used for quite a long time please switch it off and remove the plug from the electric source.

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Voltage: 220V / 110V / 50Hz. Power: 15W

Temperature: 45°C – 85°C ± 4°C.

Motor speed: 9000 rpm ± 20%

【 Important points to notice 】

In order to ensure using the Yu Yang Cup properly, please do read the below instructions carefully first.

1.Keep away from children.

2.Avoid being broken: The scraping head of Yu Yang Cup is made of ceramic. Do not throw or drop down on hard substances. Stop using the Cup if the ceramic is broken to prevent from cutting the skin.

3.Keep away from water: Clean with soft brush and antimicrobial product specially made by our Company or use a brush and cotton cloth moistened with sterile water to do the cleaning.

4.Always place the Cup on the wooden base provided

before and after the treatment to avoid getting scalded.

5. Turn off and disconnect the electricity supply when not in use.

6. Stop using the Cup if the Cup is damaged, broken or not functioning. Please contact the suppliers immediately for repairing. The Cup must be repaired by our technicians and the Company does not hold any responsibilities if the Cup is not repaired or opened by our Staff.

7. Do not use the Cup in heavy moisture environment or easily explosive surroundings.

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
1. For those who have skin festers, serious heart diseases, high blood pressure or patients with heart pacemakers, thromboembolism, purpura, leukemia, and hemophilia should avoid using the Cup or may consult their doctors before use.

2. Acute contusions such as sports injuries like sprain or bruise, should be treated with ice compress within 2 hours. Do not apply any hot compress when swelling or bleeding occurs. Otherwise the injuries would become worse. The hot compress can be used after swelling disappears.

3. Use with cautions: Pregnant women, female in menstruation period, and children should get guidance from their doctors before use. Do not apply the Cup on patients who are too exhausted, hungry, or drunken.

【 The Principle of Guasha 】

We sincerely introduce a safe and environmental Chinese therapy to you. It is called “Guasha” which means “scraping the skin”. Guasha is proved to be one of the effective but without any side effect treatments in Chinese Traditional Medicine. During the treatment some oil is applied to the part of the skin where you feel uncomfortable. Then the therapist starts “scraping” the area with a “scraper”. The scraper is made of jade or stone. Modern method is using Yu Yang Cup to scrape the skin. Yu Yang Cup can produce heat from temperature 45°C to 85°C. The warmth makes guasha much more effective. Very soon the skin becomes red. Red spots, toxic substances (痧), appear in the skin (not on the skin). The skin is very safe but no hemorrhage happens. Because the scraping is on the oily skin the patient will not feel pain – neither during nor after the treatment. The red spots will entirely disappear in 3 to 7 days.

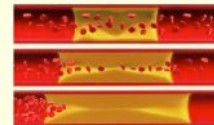Why can guasha remove toxic substances from your body?

The toxic substances including cholesterol, calcified compound, dead cells etc. inside your body accumulates on the wall of the capillaries and results in blood stagnant. It means that there are difficulties in the movement of blood.

By means of guasha, toxic things can be pressed out of the capillaries and go to the lymphatic system. The lymph fluids will then digest, engulf and dissolve these oily toxic substances into soluble substances. These soluble substances can then be excreted out of the body during urination.


【痧象圖解 Relation of red dots with degree of illness 】

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
Very little red dots show
that the body is healthy.

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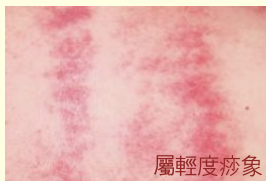Lightly dense red dots show
that the body is a bit unhealthy

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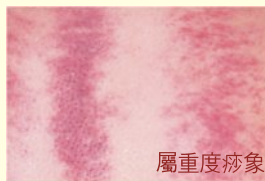
Large area of dense purple dots
exceeding 2cm diameters show
that the meridians are seriously
blocked, the patient is weak. The
illness is serious and the period
of recovery is long.

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Dense dim red dots show that the
patient is very weak and suffering
from chronic illnesses.

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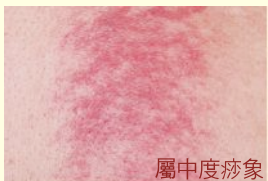
Red dense red dots show that
the body is suffering from a
temperature and inflammation.

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
Dense red dots exceeding 2cm
diameters show that the
meridians are lightly blocked.
The illness is not serious and
the period of recovery is short.

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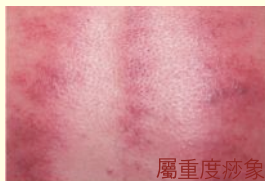Large area of dense dim purple
and red dots shows that the blood
is in stagnant condition. The
illness is serious and the period of
recovery is long.

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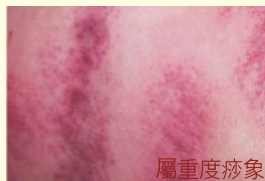Large area of dense dim purple
and dark red dots show that the
blood is seriously stagnant. The
illness is very serious and the
period of recovery is long.

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Dear customer,

Thank you for choosing the product of the company of Fu Yang
 Tong. This is a maintenance card for quality security. It is free of charge
 for repairing any malfunctioning within one year from the date of
 purchasing . For those exceeding the warranty period or the Cup is
 damaged by improper operation, the cost of repairing fee and the delivery
 fee will be charged.

Please go to your local distributor with the receipt, the warranty card and
 the Cup for repairing.

保修卡 Repair card

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維修記錄/Repair Registration

| 報修日期 Repair date | 修理日期 Repair date | 故障現象 Defect description | 故障原因 Damaged by | 修理詳情 Fix content | 修理員 Technician |
|---------------------|---------------------|-------------------------------|--------------------|---------------------|-------------------|
| | | | | | |
| | | | | | |
| | | | | | |

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Certificate

產品名稱: 御陽罐/Yu Yang Cup
Product

產品型號: YYG 001
Serial Number

檢驗員: 合格/Pass
001
Inspector

生產日期: _____
Manufacturing Date