碳纖維發熱,不含電磁輻射! Carbon fibre produces heat thus no magnetic radiation.

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Using low electrical voltage insures much safer in using our Products.

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MANUFACTURER:GUANGZHOU YUYANG TECHNOLOGI LIMITED

地 址: 廣州市番禺區南村鎮江南工業一區二横路 12 號雅景創意園 B 棟 502 室
Room 502, Building B, Yajing Creative Park, No.12, 2nd Road, Jiangnan Industrial
Zone, Xingye Avenue, Nancun, Panyu District, Guangzhou

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Your Health Partner for Better Health and Better Life

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
Body temperature affects your health and blood circulation is influenced by the coldness and warmth of the body.

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A.Characteristics of carbon fibre: Carbon fiber is inflammable and unbendable, save energy and electricity, radiate infra–red light of wave length $8{\text -}14~\mu\text{m}$ which is similar to the wave length of solar light. This light can help relieve the pain of rheumatism and periarthritis.

B. 砭石特點: 天然砭石可安神、調理氣血、疏通經絡作用。加熱後的砭石發出許多有益遠紅外線和超聲波脈衝,促進血液微循環,調理新陳代謝。B.Characteristics of Bianstone: Natural Bianstones can ease one's mind, regulate blood and Qi energy and help circulation of meridians. When Bianstones are heated, they can produce more infra-red light and ultrasonic pulses that can help the circulations of the blood capillaries to regulate the body metabolism.

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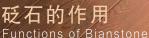Heat Treatment: improve our cell metabolism, speed our blood movement and reduce stagnant plague in our blood vessels as well as eliminate coldness and wet. 2. 遠紅外線: 有奇異的能量場, 遠紅 外線波長 8 ~ 14 微米, 可在人體產生 遠紅外線熱像並循環經絡, 增強組織活 力, 提高免疫力。

Infra-red wave: an amazing energy field of far infra-red of wave length between 8 to 14 micrometer can benefit our circulation of meridians, increase body tissue energy and improve our immune system.

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Bianstone possesses thirty or more rare elements such as Strontium, Titanium,. Chromium, Manganese and Zinc. When these elements contact our skin, they can activate the blood circulation inside our blood capillaries.

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Yu Yang Bianstone sheet heating pad is good for all four seasons. During spring and summer, your body heat initiated the energy from the stone. Plug into power outlets during autumn and winter time can also be used as hot pad to warm your body in addition to improve your health and ease all pains.

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Warming Du Meridian: Blood runs more smoothly from the waist to lower spine can sooth backache. Warming the bladder from Huiyang(BL11) can balance all internal organs in human body.

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Warming Ren Meridian: From Qugu(RN2) to Tiantu(RN22) may adjust ladies's blood flow and menstruation. Improving chance of pregnancy and skin texture.

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Warming Legs: Soothing pains, numbness of the legs. Correcting fascia to normal conditions.

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
Yu Yang Bianstone sheet heating pad: The above heating pad uses carbon fibre to produce hotness with no radiation. The stone energy produces infrared rays that can activate meridians and blood circulations thus eliminates toxic substance inside the human body for better health and relieve backaches.

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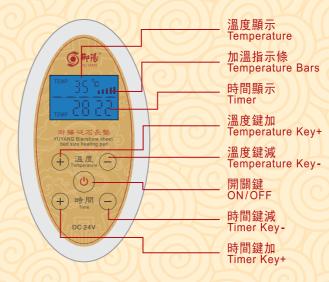1. Children and Senior citizens must be supervised by adults when use. Temperature should be adjusted to 50 degrees or below to avoid scalding.

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2. This product should never be immersed or washed in water.

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3. This product should not be used in high humidity surroundings.

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4.This product should not be used after drinking wine.

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5.Pregnant women or during menstruation period are advised not to use this product.

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6.When this product is out of order please hand it to our Company for repdir. For your safety this product should never be repaired by non-professionals.

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7.For those people who are suffering from heart aches or irregular heart pulses, please consult your doctors or medical practitioners before use.

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8.The temperature of this product ranges from 40 ~70 degree Celsius. Users can adjust the temperature to suit their own feelings.

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Reminder: The temperature shown on the pad is the temperature of the carbon fiber, the surface temperature of the pad may be affected by the environment and may have 5~10 degree difference. This difference is normal and acceptable!

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To turn on , press " o " for one second. Wait for 10 minutes for default temperature to reach 50 °C. Difference will be + or – 5 °C. Default time set for 60 minutes will then countdown to zero. Press " o " for one second to turn off. Unplug when not in use.

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Temperature setting: Press temperature button " \oplus " or " \ominus " for desired temperature. Desired temperature will be set after releasing the button after three seconds. Temperatures are set as 40/45/50/55/60/65/70 degrees.

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Time setting: Press Time button " ⊕ " or " ⊝ " for desired time. Time set will countdown to zero. Time is set as 60/90/120 minutes

4. 開機後,顯示屏工作中顯示的溫度是實時溫度,不是顯示設定溫度。 The temperature shown on the panel when the power is on is the room temperature. The tool will reach the desired temperature after a few minutes

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
Sound indications: All buttons will beep once when pressed. You will hear one short beep after temperature or time is set and five beeps when time ends and turning off automatically. Long nonstop beeping means it is not functioning well.

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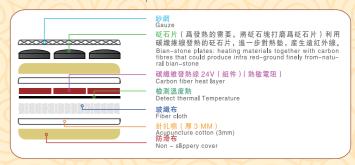When the display panel shows E1, E2, E3, please contact our company immediately for repairs.

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